English: Triplets Cable for Series 50 Fetal Monitors
Instructions for Use

- 1 x M1360-61671 triplets cable
- 1 x Fetal Monitoring Documentation CD-ROM
- These instructions

The M1360-61671 triplets cable lets you connect either an 8040A or an 8041A Fetal Monitor to a M1350A (Series 50 IX) or a M1350B/C (Series 50 XM/XMO) Intrapartum Fetal/Maternal Monitor. This connection lets you monitor three non-invasive fetal heart rates (FHRs) plus one maternal heart rate (MHR).

To use this connection you need:


- M1350A or M1350B/C with telemetry interface or telemetry system interface.
- Standard 8041A or 8040A with interface option J01.

The M1350A or M1350B/C measures the FHR of triplets 1 and 2. The 8040A/8041A measures the FHR of the third triplet and sends the information to the M1350A/B/C ("master device"). All three measurements are documented on the M1350A/B/C.

The third trace is not transmitted to OB TraceVue or any other obstetrical surveillance system, so it is important to retain the paper trace as part of the medical records.


Installation should be carried out by qualified technical personnel.

Connect the cable provided to the 25 pin connector on the telemetry interface or telemetry system interface of the M1350A or M1350B/C as shown in the following figure:

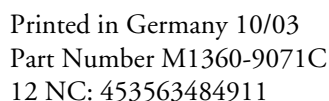
Connect the other end of the cable to the system interface of either the 8040A or the 8041A as shown in the following figures:

8040A
rear view

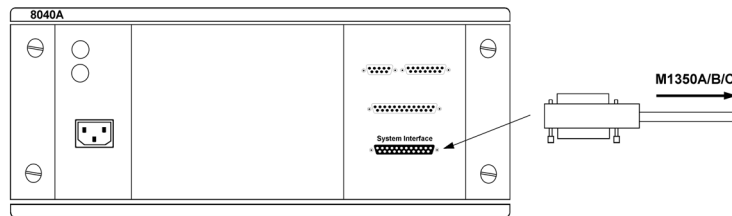8041A
rear view

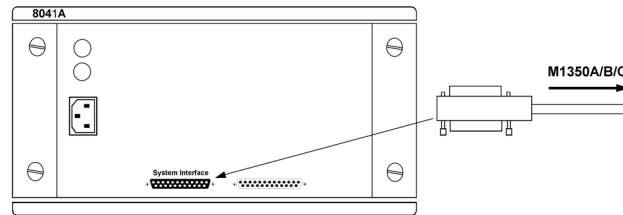The equipment should already have been tested according to the safety tests outlined in the M1350B/C *Instructions for Use* and detailed in the M1350B/C *Service Guide*, both contained on the CD-ROM that comes with the cable.

After connecting the equipment together, perform the system test, if applicable, and the performance test for the M1350A/B/C (see the M1350B/C *Instructions for Use* on the CD-ROM that comes with the cable. This is also available in printed form on demand. Please contact your local Philips representative).

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